




## Multi-sample soil respiration system

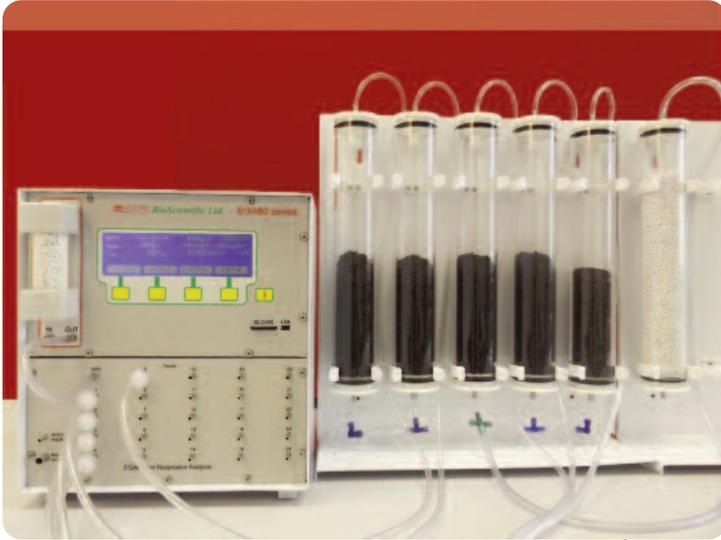

- + Dual CO<sub>2</sub>/H<sub>2</sub>O analysis
- + Up to 24 soil samples with one unit
- + Accurate and proven technology
- + Compact and fully integrated
- + Flow maintained to all selected channels
- + Flexible channel selection
- + SD card data storage and USB output

### ADC BioScientific Ltd.

World leaders in developing and manufacturing multi-sample analysis systems (such as the renowned ADC225), where a number of soil samples may be multiplexed to and measured by a single CO<sub>2</sub> analyser.

### Compact and fully integrated

ADC now introduces the EGA60, a fully integrated, soil respiration system, featuring an accurate and reliable CO<sub>2</sub> analyser, combined with a versatile gas multiplexer in one compact unit. The EGA60 is designed for long-term, continuous experiments.

### Analysis of up to 24 soil samples

A single EGA60 system can sequentially analyse up to 24 soil samples. The EGA60 is available with 5, 10, 15, 20 or 25 channels. One channel being reserved for an autocalibration 'zero' column.

### Fully programmable

Complete functionality is achieved with just 5 keys driving a series of menus. Sampling times of individual channels can be set, together with total experimental times. Settings can be copied over to several channels, saving time.

### Constant flow maintains sample integrity

A constant flow is maintained around the system to all sample channels, at all times. This prevents the build up of high CO<sub>2</sub> concentrations within any sample chamber. Flow in each channel may be automatically programmed.

### Versatile

EGA60 systems are available with 500ml cylindrical, glass or acrylic, soil columns. Packages are available with 5, 10, 15, 20 or 25 columns attached to adaptable mounting plates.

### Applications:

Soil biomass determination (volume of CO<sub>2</sub> released per unit mass of soil over time).

Toxicology studies, agricultural registration (ISO 14240-1).

Plant physiology, environmental chambers, insect respiration and fruit storage

## Integral data storage

The EGA60 provides integral data storage on interchangeable SD cards, each capable of storing many hundreds of thousands of data points.

Data may be downloaded to a PC directly from the SD card or via a USB port. Data opens in a comma delimited format as an Excel spreadsheet.

## Determination of soil biomass

The EGA60 can also be configured with your own columns or chambers, whether your application is soil ecology or soil treatment.

Data provided by the EGA60 can subsequently be used to calculate the volume of CO<sub>2</sub> released per unit mass of soil over time.

## Soil toxicology

One of the main applications for the EGA60 is in agrochemical registration:-

Prior to agrochemical registration for commercial usage vigorous toxicology tests must be carried out. This includes active, aerobic, heterotrophic microbial biomass degradation in aerated agricultural and mineral soils. To perform these tests an international standard has been published.

The ISO 14240-1 Soil quality - Determination of soil microbial biomass, Part 1: Substrate-induced respiration method (SIR) was based on measurements made using an earlier ADC soil respiration system (Heinemeyer *et al.*, 1989).

## Technical Specification

### Measurement range and technique:

CO<sub>2</sub>: 0-2000ppm, 1ppm resolution Infrared gas analysis

H<sub>2</sub>O: 0-75mbar, 0.1mbar resolution.

Two laser - trimmed, fast response water vapour sensors

**Flow control:** 0 to 500ml min<sup>-1</sup> on each channel

**Test duration:** Set by time or number of measurement cycles

**Dwell time:** 2 seconds to 999 minutes on each channel

**Warm up time:** 5 minutes @ 20°C

**Display:** 240 x 64 graphic LED backlit LCD

**Recorded data:** Removable SD card (1Gb card typically stores 16 million sets of data)

**Power supply:** 230/110V 50/60Hz

### Electrical outputs:

USB connection: Mini-B  
RS232: 9 Pin "D" type

**Analogue output:** 0-5V or 4-20mA

**Analogue inputs:** Seven 0-5V

**Operating temperature range:** 5°C to 45°C

**Dimensions:** 27 x 25 x 15cm

**Weight:** 7.5kg

## ADC Never Compromise on Quality

*"Quality of product and quality of service."*

From design to delivery, ensuring optimal performance and reliability is of paramount importance to our team of experienced engineers. Once in the field, you are supported by our worldwide network of customer support centres.



Büro Cottbus

Taubenstrasse 4  
D - 03046 Cottbus  
vertrieb@upgmbh.com

Fon: +49 355 485 54-0  
Fax: +49 355 485 54-15  
www.upgmbh.com

ADC Bioscientific Ltd., Global House, Geddings Road,  
Hoddesdon, Herts EN11 0NT, UK

Web: [www.adc.co.uk](http://www.adc.co.uk) Email: [sales@adc.co.uk](mailto:sales@adc.co.uk) Tel: +44 (0)1992 464527